


# South West London Wound Management Products Formulary - January 2019

A collaboration between acute trusts, clinical commissioning groups & community health providers in south west London

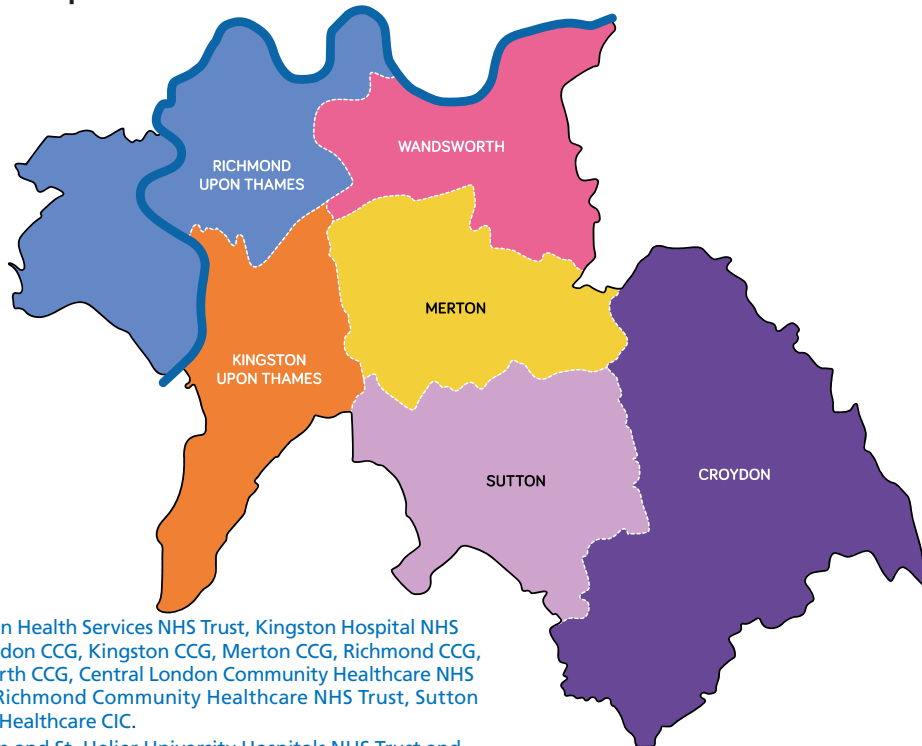

**Approved by:** Croydon Health Services NHS Trust, Kingston Hospital NHS Foundation Trust, Croydon CCG, Kingston CCG, Merton CCG, Richmond CCG, Sutton CCG, Wandsworth CCG, Central London Community Healthcare NHS Trust, Hounslow and Richmond Community Healthcare NHS Trust, Sutton Health and Care, Your Healthcare CIC.

**Supported by:** Epsom and St. Helier University Hospitals NHS Trust and St. Georges University Hospitals NHS Foundation Trust for use by primary care and community health providers

## Introduction to the formulary:

- This formulary has been developed in response to an identified need across the health care sector in south west London. It includes recommended dressings and products for the management of wounds in south west London. The choice of product is based upon national guidance, research based evidence, practical experience and cost, and has been recommended by the South West London Wound Care Sub-Group which includes tissue viability specialists from each organisation.
  - The formulary will be updated periodically by the wound care sub-group should any changes be required to reflect best practice (e.g. an update to national guidelines). Please contact your local tissue viability nurse if you would like to discuss any possible changes to the formulary. The South West London Wound Care Product Evaluation Standard Operating Procedure will be followed for the addition of new products onto the formulary. The most up to date version of the formulary and other wound care resources will be available electronically on the south west London website available here <http://www.swlmcg.nhs.uk/Clinical/Pages/Wound-Care.aspx>
  - The aim of this formulary is to guide the clinician in their choice of dressings so that an optimum wound healing environment can be achieved whilst using cost effective treatment options, with minimal disruption to the patient between care settings. It should also:
    - Encourage good clinical practice that is research and evidence based
    - Provide healthcare professionals with an up to date reference tool
    - Promote a rational approach towards prescribing wound care products across the healthcare sectors
    - Reduce confusion across healthcare settings
    - Reduce wastage
  - This formulary is to be used in conjunction with local policies and holistic assessment, and practitioners should use their own professional judgement in individual patient circumstances. This is a stand-alone formulary to guide choice of products – for wound assessment see separate local wound management guidelines.
  - The clinician should read all product literature and manufacturer's guidelines prior to use. The advice provided in this formulary is not comprehensive and it is the prescriber's responsibility to consult other appropriate sources where necessary.
  - This formulary does not include all items that are available on prescription.
  - Specialist products highlighted in red must only be used following consultation with a specialist in wound care.
  - Dressings may be prescribed as single items and not just by pack/box. Please prescribe the minimum number of dressings required as wound healing is a dynamic process and the characteristics of a dressing required by wounds and wound size can change over time.
  - The tissue viability team for your area is available to discuss dressing selection. Contact details are at the end of this formulary.
-

# Contents

# Page

Adhesive dressings for lightly exuding wounds	2
Alginate dressings	2
Antimicrobial dressings – honey	2
Antimicrobial dressings – iodine	3
Antimicrobial dressings – silver	3
Antimicrobial dressings – other	4
Atraumatic and low adherence dressings	4
Bandages	5, 6, 7
Foam dressings	8, 9
Hydrocolloid dressings	10
Hydrofibre dressings	10
Hydrogel dressings	11
Negative pressure wound therapy (npwt)	11
Odour absorbent dressings	12
Stockinettes	12
Super absorbent dressings for moderate to heavily exuding wounds	13
Tapes	13
Vapour permeable film dressings	14
Miscellaneous	15, 16
Contact details for local tissue viability service	17

PRODUCT	SIZE	COMMUNITY COST PER DRESSING*	COMMENT
---------	------	------------------------------	---------

### ADHESIVE DRESSINGS FOR LIGHTLY EXUDING WOUNDS

Absorbent dressings with adhesive border for lightly exuding wounds. Can be used as primary or secondary dressings and have limited absorbency.

Softpore	6cm x 7cm	£0.06	
	10cm x 10cm	£0.12	
	10cm x 20cm	£0.35	
	10cm x 30cm	£0.49	

### ALGINATE DRESSINGS

Has haemostatic properties. Only use on bleeding wounds. These dressings are NOT suitable for use on dry wounds. Cavity wounds should not be over-packed as the fibres swell and can lead to tissue ischaemia. Cut sheets to the size of the wound to avoid maceration and excoriation of the surrounding skin, as any exudate will transfer laterally.

ActivHeal Alginate	5cm x 5cm	£0.59	Irrigate with saline or water to remove if the dressing is adhered to the wound. Cut sheets to size of the wound to avoid maceration of surrounding skin.
	10cm x 10cm	£1.15	
	10cm x 20cm	£2.83	

### ANTIMICROBIAL DRESSINGS – HONEY

Suitable for critically colonized or locally infected wounds. Systemic antibiotics are indicated in cases where there are clinical signs of systemic infection although in patients with diabetes these symptoms could be masked. The tissue viability team should be contacted if use of antimicrobials is considered for more than 14 days.

Medihoney Antibacterial Honey Apinate	5cm x 5cm	£2.01	Effective against a broad spectrum of bacteria. Fast, effective autolytic debridement on sloughy & necrotic tissue, removes malodour, provides moist wound environment. Less complaints of stinging/drawing effect. Contraindications - patients with known allergies to bee venom.
	10cm x 10cm	£3.42	
Medihoney Antibacterial Honey Apinate rope	1.9cm x 30cm	£4.23	
Medihoney HCS	6cm x 6cm	£2.26	
	11cm x 11cm	£4.50	
Medihoney HCS Surgical with Adhesive Border	7.5cm x 20cm	£3.07	
Medihoney HCS with Adhesive Border	11cm x 11cm	£3.08	
	15cm x 15cm	£5.83	
Medihoney Tulle dressing	5cm x 5cm	£1.71	
	10cm x 10cm	£3.00	
Medihoney Antibacterial wound gel	10g	£2.71	As above. Reusable for the same patient for up to 4 months.
	20g	£4.05	

PRODUCT	SIZE	COMMUNITY COST PER DRESSING*	COMMENT
---------	------	------------------------------	---------

#### ANTIMICROBIAL DRESSINGS – IODINE

Suitable for critically colonised or locally infected wounds. Systemic antibiotics are indicated in cases of overt wound infection where the classical signs are evident. The tissue viability team should be contacted if use of antimicrobials is considered for more than 14 days.

Inadine	5cm x 5cm	£0.28	For lightly exuding wounds. Contraindications: patients receiving lithium, in thyroid disorders, in pregnancy and breast feeding, and in children; they should be used with caution in patients with severe renal impairment or history of thyroid disorder.
	9.5cm x 9.5cm	£0.42	
Iodoflex paste	5g	£4.19	Do not use on patients with sensitivity to any of products ingredients, children, pregnant or lactating women or those with thyroid disorders or renal impairment
	10g	£8.37	
	17g	£13.25	

#### ANTIMICROBIAL DRESSINGS – SILVER

Suitable for critically colonized or locally infected wounds. Systemic antibiotics are indicated in cases of overt wound infection where the classical signs are evident. Not to be used on patients with a known sensitivity to Silver. The tissue viability team should be contacted if use of antimicrobials is considered in secondary care or for more than 14 days in primary care.

Aquacel Ag+ Extra	5cm x 5cm	£1.98	Review with TVN if required for more than 14 days in primary care. - please do not cut on the stitches line - at least 2cm overlap from the peri-wound
	10cm x 10cm	£4.71	
	15cm x 15cm	£8.88	
	4cm x 10cm	£2.87	
Aquacel Ag+ ribbon	1cm x 45cm ribbon	£3.10	
	2cm x 45cm ribbon	£4.74	
Acticoat Flex 3	5cm x 5cm	£3.60	For use by specialists in wound care only or under their advice and direction. Dressing can stay in place for 3 days.
	10cm x 10cm	£8.80	
	10cm x 20cm	£13.75	
Actisorb silver 220	6.5cm x 9.5cm	£1.66	For use by specialists in wound care only or under their advice and direction.
	10.5cm x 10.5cm	£2.60	

PRODUCT	SIZE	COMMUNITY COST PER DRESSING*	COMMENT
---------	------	------------------------------	---------

#### ANTIMICROBIAL DRESSINGS – OTHER

Suitable for critically colonized or locally infected wounds. Systemic antibiotics are indicated in cases of overt wound infection where the classical signs are evident. The tissue viability team should be contacted if use of antimicrobials is considered for more than 14 days.

Cutimed sorbact swab	4cm x 6cm (11cm x 6cm unfolded)	£1.68	Binds bacteria to the swab. Requires a secondary dressing.
	7cm x 9cm (17cm x 20cm unfolded)	£2.80	
Flaminal Forte	15g tube (approx. coverage 40cm <sup>2</sup> )	£7.79	<b>For use by specialists in wound care only</b> or under their advice and direction.  For wounds that have moderate to heavy amounts of exudate. Can be used throughout all the stages of wound healing. A protease modulator.
	50g tube (approx. coverage 130cm <sup>2</sup> )	£25.78	
Flaminal hydro	15g tube (approx. coverage 40cm <sup>2</sup> )	£7.79	<b>For use by specialists in wound care only</b> or under their advice and direction.  For wounds that have mild to moderate amounts of exudate. Can be used throughout all the stages of wound healing. A protease modulator.
	50g tube (approx. coverage 130cm <sup>2</sup> )	£25.78	

#### ATRAUMATIC AND LOW ADHERENT DRESSINGS

Used to prevent trauma to granulating or friable wounds.

Profore contact layer WCL	14cm x 20cm	£0.33	Cut to size
Atrauman	5cm x 5cm	£0.34	Effective for up to 7 days
	7.5cm x 10cm	£0.35	
	10cm x 20cm	£0.80	
	20cm x 30cm	£2.20	
Silflex	5cm x 7cm	£1.12	Indications include skin tear, burns and surgical wounds Maximum wear time of 14 days.
	8cm x 10cm	£2.29	
	12cm x 15cm	£4.61	
NA Ultra	9.5cm x 9.5cm	£0.34	
	19cm x 9.5cm	£0.64	

PRODUCT	SIZE	COMMUNITY COST PER DRESSING*	QTY PER PACK	COMMENT
---------	------	------------------------------	--------------	---------

### BANDAGES AND COMPRESSION SYSTEMS

Compression bandages should only be applied following a holistic leg ulcer assessment inclusive of Doppler ultrasound. Compression bandages should be selected carefully according to ankle circumference and limb length. Compression bandages are contraindicated in arterial leg ulcers. **Compression bandages should only be applied by staff that have completed the local leg ulcer management training programme and they must be signed off as competent in the application of compression bandages. Wraps can be applied by anyone that has been shown how to use as part of self-care e.g. carer, patient.** Only order quantities needed.

Actico	8cm x 6m	£3.31	1	Cohesive inelastic short stretch bandage, requires undercast padding to protect bony prominences.
	10cm x 6m	£3.43	1	
	12cm x 6m	£4.38	1	
Peha-haft – latex free; double adhesive effect	2.5cm x 4m	£0.75	1	
	4cm x 4m	£0.48	1	
	6cm x 4m	£0.57	1	
	8cm x 4m	£0.68	1	
	10cm x 4m	£0.78	1	
	12cm x 4m	£0.92	1	
Multi-layer compression bandage kit, Urgo K-four:				Prescribe the individual bandages below, not the kit, to reduce waste.
K-soft	10cm x 3.5m	£0.46	1	Soft absorbent bandage for redistribution pressure over bony prominences. Applied spiral with 50% overlap.
K-soft long	10cm x 4.5m	£0.58	1	
K-lite	10cm x 4.5m	£1.02	1	Light Support type 2 bandage. Applied spiral with 50% overlap.
	15cm x 4.5m	£1.47	1	
K-lite long	10cm x 5.25m	£1.16	1	Light compression type 3a bandage. Applied in figure of 8 with 50% stretch.
K-plus	10cm x 8.7m	£2.31	1	
K-plus long	10cm x 10.25m	£2.67	1	
KO-Flex	10cm x 6m (stretched)	£3.07	1	Cohesive bandage. Applied spiral with 50% overlap and 50% stretch
KO-Flex long	10cm x 7m (stretched)	£3.51	1	
K-ThreeC	10cm x 3cm (unstretched)	£2.87	1	

CONTINUES ON NEXT PAGE

PRODUCT	SIZE	COMMUNITY COST PER DRESSING*	QTY PER PACK	COMMENT
UrgoKTwo 2 layer compression bandaging kit.	Ankle circumference			To ensure the correct K-Two kit is prescribed it is necessary to write the required bandage width after the size of kit as 8cm, 10cm and 12cm widths available.
	18-25cm (8cm)	£7.77	1	
	18-25cm (10cm)	£8.24	1	
	18-25cm (12cm)	£10.39	1	
	25-32cm (8cm)	£8.45	1	
	25-32cm (10cm)	£9.00	1	
UrgoKTwo Reduced - multi-layer compression bandage kit	18-25cm	£8.24	1	
	25cm - 32cm	£9.01	1	
<p>Comprising: K-Tech – soft absorbent bandage for redistribution pressure over bony prominences. Applied spiral with 50% overlap for 18cm-25cm and 2/3 overlap for 25-32cm and 50% stretch.</p> <p>K-Press – Cohesive bandage. Applied spiral with 50% overlap for 18cm- 25cm and 2/3 overlap for 25-32cm and 50% stretch.</p>				
Juxta cures compression wraps	Short, Standard, Long	£156.20	1 to last up to 6 months	Juxtacures provides below the knee compression and is used in the management of venous leg ulcers (VLUs) when oedema is present. Guaranteed for up to 6 months for leg ulcer management.
Juxta fit compression wraps	Lower Legging 28cm, S, M, L, XL, XXL	£136.90	1 to last up to 12 months	Juxtafit products are used for lymphoedema, and come in a variety of styles to provide arm or leg compression as required.  Guaranteed for up to 12 months for chronic oedema/ lymphedema management.
	Lower Legging 36cm S, M, L, XL, XXL	£136.90		
	Lower Legging Extra Wide Calf 28cm M, L, XL	£132.66		
	Lower Legging Extra Wide Calf 36cm M, L, XL	£132.66		
Ankle Foot Wrap, Interlocking	S, M, L	£38.33		
Ankle Foot Wrap, Closed Heel	S, M, L	£38.33		

CONTINUES ON NEXT PAGE



PRODUCT	SIZE	COMMUNITY COST PER DRESSING*	QTY PER PACK	COMMENT
Juxta lite compression wraps	Available as Short and Long.  S, M, L, XL, XXL  Available as Short Extra Wide and Long Extra Wide.  M, L, XL,	£97.72	1 to last up to 6 months	Juxtalite also provides below the knee compression, and is suitable where there is no oedema and for the prevention of recurrence.  Guaranteed for up to 6 months for leg ulcer prevention or maintenance post healing.  They can be applied by the patient themselves, carers and clinicians. Better patient concordance with compression.
Ichthopaste	7.5cm x 6 m	£3.78	6	Zinc paste and Ichthammol bandage. Patch test for skin sensitivity for 4 days before use - please refer to packet insert for application technique.
Viscopaste PB7	7.5cm x 6 m	£3.74	6	Zinc paste bandage. Patch test for skin sensitivity for 4 days before use -please refer to packet insert for application technique.
Mollelast bandage	4cm x 4m	£0.31	1	<b>For use by specialists in wound care only</b> or under their advice and direction. For the management of toe oedema.

PRODUCT	SIZE	COMMUNITY COST PER DRESSING*	COMMENT
---------	------	------------------------------	---------

#### NON-ADHESIVE FOAM DRESSINGS

Change according to clinical need. Creams and ointments can reduce the adherence of some products. Do NOT cover with film dressings.

<b>Advazorb</b>	5cm x 5cm	£0.65	Primary or secondary dressing for low to moderately exuding wounds. Non-adhesive.
	7.5cm x 7.5cm	£0.79	
	10cm x 10cm	£1.09	
	10cm x 20cm	£3.37	
	12.5cm x 12.5cm	£1.60	
	15cm x 15cm	£2.12	
	20cm x 20cm	£3.78	
<b>Advazorb heel</b>	17cm x 21cm	£4.78	
<b>Urgotul Absorb</b>	6cm x 6cm	£1.21	Primary or secondary dressing for moderately exuding wounds. Non-adhesive.
	10cm x 10cm	£2.38	
	15cm x 20cm	£4.20	
<b>Biatain Non-Adhesive</b>	5cm x 7cm	£1.34	Primary or secondary dressing for moderately exuding wounds. Non-adhesive.
	10cm x 10cm	£2.44	
	15cm x 15cm	£4.50	
	10cm x 20cm	£4.03	

#### ADHESIVE FOAM DRESSINGS

Change according to clinical need. Creams and ointments can reduce the adherence of some products. Do NOT cover with film dressings.

<b>FoamLite</b>	5cm x 5cm	£0.75	Primary or secondary dressing for low to moderately exuding wounds. Adhesive.
	8cm x 8cm	£1.15	
	10cm x 10cm	£1.57	
	10cm x 20cm	£2.05	
	15cm x 15cm	£3.50	
	5.5cm x 12cm	£1.56	

CONTINUES ON NEXT PAGE

PRODUCT	SIZE	COMMUNITY COST PER DRESSING*	COMMENT
Urgotul Absorb Border	6.5cm x 10cm	£1.48	Primary or secondary dressing for moderately exuding wounds. Adhesive.
	8cm x 8cm	£1.38	
	8cm x 15cm	£2.44	
	10cm x 10cm	£1.98	
	10cm x 25cm	£3.67	
	13cm x 13cm	£2.45	
	15cm x 20cm	£3.98	
Urgotul Absorb Border Sacrum	20cm x 20cm	£4.35	
Biatain Silicone	7.5cm x 7.5cm	£1.48	<b>Second Line for fragile skin</b> Primary or secondary dressing for low to moderately exuding wounds with sensitive/fragile surrounding skin.
	10cm x 10cm	£2.18	
	10cm x 20cm	£2.75	
	10cm x 30cm	£5.29	
	12.5cm x 12.5cm	£2.67	
	15cm x 15cm	£3.96	
Biatain Silicone Multishape	14cm x 19.5cm	£4.39	

#### FOAM DRESSINGS FOR WOUNDS ALSO REQUIRING A HYDROFIBRE DRESSING

Change according to clinical need. Creams and ointments can reduce the adherence of some products. Do NOT cover with film dressings.

Aquacel Foam Non-adhesive	5cm x 5cm	£1.37	<b>ONLY use when a hydrofibre dressing AND foam dressing would otherwise be used.</b> Primary or secondary dressing for moderate to high exuding wounds. Absorbent foam dressing with integral Aquacel wound contact layer.
	10cm x 10cm	£2.60	
	15cm x 15cm	£4.37	
	10cm x 20cm	£3.59	
Aquacel Foam Adhesive	8cm x 8cm	£1.42	
	10cm x 10cm	£2.20	
	12.5cm x 12.5cm	£2.72	
	17.5cm x 17.5cm	£5.45	
Aquacel Foam Adhesive (Heel)	19.8cm x 14cm	£5.57	

PRODUCT	SIZE	COMMUNITY COST PER DRESSING*	COMMENT
---------	------	------------------------------	---------

### HYDROCOLLOID DRESSINGS

Self-adhesive and waterproof dressings for the management of light to moderately exuding wounds. Facilitate rehydration and autolytic debridement of dry, sloughy or necrotic wounds. Also suitable for promoting granulation tissue. Hydrocolloid wafer dressings are not recommended for heavily exuding or infected wounds. Hydrocolloid wafer dressings should extend at least 2cm beyond the edge of the wound. **Hydrocolloids are not recommended for use on diabetic foot wounds as they are too occlusive.**

Duoderm Extra Thin	7.5cm x 7.5cm	£0.82	For dry to lightly exuding wounds. Not for use in diabetic foot ulcer, otherwise seek TVN advice. Please consider using adhesive remover for fragile skin.
	10cm x 10cm	£1.36	
	15cm x 15cm	£2.93	
Duoderm signal	10cm x 10cm	£2.17	Has a border. For dry to lightly exuding wounds. Not for use in diabetic foot ulcer, otherwise seek TVN advice. Please consider using adhesive remover for fragile skin.
	14cm x 14cm	£3.81	
Gelatine-free hydrocolloid:			
Tegaderm hydrocolloid thin	10cm x 12cm	£1.56	For dry to lightly exuding wounds. Not for use in diabetic foot ulcer. Please consider using adhesive remover for fragile skin.
	13cm x 15cm	£2.92	
Tegaderm hydrocolloid	10cm x 10cm	£2.39	

### HYDROFIBRE DRESSINGS

Hydrofibre is a fibrous hydrocolloid. It has the ability to absorb and lock away large amounts of exudate; it prevents lateral wicking of exudate and forms a soft gel. Not suitable for use on dry wounds. Do not use in combination with hydrogels. Requires a secondary dressing.

KerraCel	5cm x 5cm	£0.60	Should overlap 1cm onto the skin surrounding the wound.
	10cm x 10cm	£1.44	
	15cm x 15cm	£2.71	
Aquacel ribbon	1cm x 45cm	£1.89	For cavity/sinus wounds.
Aquacel ribbon	2cm x 45cm	£2.52	

PRODUCT	SIZE	COMMUNITY COST PER DRESSING*	COMMENT
---------	------	------------------------------	---------

### HYDROGEL DRESSINGS

A primary wound dressing product consisting of insoluble polymers which are hydrophilic and contain large amounts of water. Used primarily to promote autolytic debridement of sloughy/necrotic wounds. Can donate liquid to produce a moist environment at the surface of the wound. Not suitable for heavily exuding wounds as can lead to skin maceration/excoriation. Consider the site of the wound whether an amorphous gel (ActivHeal Hydrogel) or hydrogel sheet (Kerralite Cool) is the most suitable.

Kerralite Cool	6cm x 6cm	£1.77	Not for use on heavily exuding wounds.
	8.5cm x 12cm	£2.61	
	18cm x 12.5cm	£3.77	
KerraLite Cool Adhesive	8cm x 8cm	£2.05	
	11cm x 11cm	£2.73	
	15cm x 15cm	£4.35	
ActivHeal Hydrogel	8g	£1.23	Not for use on heavily exuding wounds. Single Use Only.
	15g	£1.41	

### NEGATIVE PRESSURE WOUND THERAPY (NPWT)

Topical negative pressure applied to the wound via foam or gauze dressing, sealed with an airtight adhesive drape which creates a vacuum via a suction device. Suitable for wounds that have large soft tissue defects, e.g. post-operative wounds, pressure ulcers, diabetic foot ulcers, leg ulcers, and fixation of skin grafts. Contraindications: Acute bleeding or risk of bleeding; exposed vessels and vascular prostheses; fistulas of unknown origin; untreated osteomyelitis; malignant wounds (except under medical guidance for palliative management of exudate); necrotic wound beds; patients with cognitive problems – review each patient individually as not suitable for all patients.

Dressings and consumables - Each area uses different brands			<p><b>For use by specialists in wound care only</b> or under their advice and direction.</p> <p>All clinicians applying NPWT must be competent in the application of NPWT.</p> <ul style="list-style-type: none"> <li>• Ensure TVN is aware of all referrals for NPWT.</li> <li>• Contact TVN if further advice/ support required</li> <li>• Follow local policy before ordering</li> </ul>
--	--	--	---

PRODUCT	SIZE	COMMUNITY COST PER DRESSING*	COMMENT
---------	------	------------------------------	---------

### ODOUR ABSORBENT DRESSINGS

Dressings containing activated charcoal that is effective in absorbing chemicals released from malodorous wounds. These can be used as primary or secondary dressings in the management of exuding, purulent and contaminated wounds (e.g. fungating) complicated by bacterial infection and offensive odour. May adhere to wounds that have dried out in which case a non-adhesive contact layer will be required. Can be used as a secondary dressing over a primary dressing to suppress odour. Moisture inactivates charcoal.

Clinisorb	10cm x 10cm	£1.94	
	10cm x 20cm	£2.58	

### STOCKINETTES

Used for dressing retention. Comfifast and Tubifast are elasticated viscose stockinettes. Consider using 5 metre packs if used regularly.

Comfifast stockinette 3.5cm (Red Line)	3.5cm x 1m	£0.56	Elasticated viscose stockinette. Use for dressing retention.
Comfifast stockinette 5cm (Green Line)	5cm x 1m	£0.58	
	5cm x 3m	£1.62	
	5cm x 5m	£2.81	
Comfifast stockinette 7.5cm (Blue Line)	7.5cm x 1m	£0.77	
	7.5cm x 3m	£2.13	
	7.5cm x 5m	£3.74	
Comfifast stockinette 10.75cm (Yellow line)	10.75cm x 1m	£1.20	
	10.75cm x 3m	£3.49	
	10.75cm x 5m	£6.04	
Comfifast stockinette 17.5cm (Beige Line)	17.5cm x 1m	£1.83	

PRODUCT	SIZE	COMMUNITY COST PER DRESSING*	COMMENT
---------	------	------------------------------	---------

### SUPERABSORBENT DRESSINGS FOR MODERATE TO HEAVILY EXUDING WOUNDS

Change according to clinical need. Use with caution on bleeding/arterial wounds as they can exacerbate the bleeding. Creams and ointments can reduce the absorbency of some of the products. Absorbency may be reduced under compression.

<b>Zetuvit Plus</b>	10cm x 10cm	£0.83	<b>First Line</b> Primary or secondary dressing for moderate to highly exuding wounds, if used under compression ensure that limb is shaped to an inverted cone shape in order that graduated compression will be achieved.
	10cm x 20cm	£1.15	
	15cm x 20cm	£1.32	
	20cm x 25cm	£1.80	
	20cm x 40cm	£2.79	
<b>KerraMaxCare</b>	5cm x 5cm	£1.01	<b>Second Line</b> Primary or secondary dressing for moderate to highly exuding wounds, thinner dressing even when used under compression, non-adhesive dressing you can use either side.
	10cm x 10cm	£1.31	
	10cm x 22cm	£1.71	
	13.5cm x 15.5cm	£1.55	
	20cm x 22cm	£3.02	
	20cm x 30cm	£3.45	

### TAPES

For fixation of dressings and bandages

<b>Hypafix</b>	5cm x 5m	£1.47	Can be applied directly onto skin.
	10cm x 5m	£2.47	
	5cm x 10m	£2.73	
	10cm x 10m	£4.76	
<b>Micropore</b>	1.25cm x 5m	£0.62	
	2.5cm x 5m	£0.92	
	5cm x 5m	£1.63	

PRODUCT	SIZE	COMMUNITY COST PER DRESSING*	COMMENT
---------	------	------------------------------	---------

#### VAPOUR PERMEABLE FILM DRESSINGS

Comfortable, waterproof film dressings that are vapour permeable allowing the passage of water vapour and oxygen. For use as a secondary dressing over hydrogels or as a primary dressing to protect fragile skin. Do not use to cover foam dressings as it reduces moisture vapour transmission rate.

Clearfilm	6cm x 7cm	£0.25	Without pad.
	10cm x 12cm	£0.55	
	12cm x 12cm	£0.62	
	15cm x 20cm	£1.15	
	20cm x 30cm	£1.95	
Clearpore	6cm x 7cm	£0.12	With pad.
	6cm x 10cm	£0.15	
	10cm x 10cm	£0.20	
	15cm x 10cm	£0.24	
	20cm x 10cm	£0.36	
	25cm x 10cm	£0.40	
Clearfilm IV	6cm x 7cm	£0.28	For use on IV/subcutaneous therapy sites.
	7cm x 9cm	£0.36	
	10cm x 12cm	£0.98	
	10cm x 15cm	£1.40	



PRODUCT	SIZE	COMMUNITY COST PER DRESSING*	QTY FOR 2 WEEK SUPPLY	COMMENT
<b>MISCELLANEOUS</b>				
Adhesive remover: Opus Lift Plus 360 Medical Adhesive Remover/Citrus Medical Adhesive remover	50ml spray	£5.99	1	Wipes should ONLY be prescribed for those that need them, e.g. for people who lack the manual strength or dexterity to use a spray.
Adapt Universal Remover Wipes	50 wipes	£14.06	1	Either a spray or wipes should be prescribed, not both.
Irrigation Solutions: Sodium Chloride 0.9% (brands change price regularly: check drug tariff for most cost-effective choice), Clinipod is currently the most cost effective option at time of production of the formulary	20ml unit dose vials	Approx. £4.80 - £7.96	15 unit dose vials	Usually 25pcs in a box.
Octenilin Solution	350ml bottle	£4.60	1 bot	<b>For use by specialists in wound care only</b> or under their advice and direction. Wound irrigation solution which aids the removal of biofilm and wound debris.
Octenillin wound gel	20ml	£4.78	n/a	<b>For use by specialists in wound care only</b> or under their advice and direction. Gel containing Betaine and PHMB that removes biofilm and wound debris. Can re-use.
For debridement: Debrisoft	10cm x 10cm	£6.55	n/a	For debridement of loose slough and debris. Single use only. Should not be left on the wound bed.
For debridement: UrgoClean Pad	6cm x 6cm	£0.97	n/a	<b>For use by specialists in wound care only or under their advice and direction.</b>  <b>Slough trapping.</b>
	10cm x 10cm	£2.16		
	20cm x 15cm	£4.06		
UrgoClean Rope (with probe)	2.5cm x 40cm	£2.43	n/a	
	5cm x 40cm	£3.21		
Footwear: Kerraped all purpose boot	S, M, L, XL	£17.92	2 max	Only for use if patient has no other suitable footwear.

CONTINUES ON NEXT PAGE

PRODUCT	SIZE	COMMUNITY COST PER DRESSING*	QTY FOR 2 WEEK SUPPLY	COMMENT
Kerrapro Pressure Reducing Pads (silicone gel pads)	Sacrum/ankle One size	£17.09	1	Polymer gel-shaped pad that redistributes pressure while protecting and padding bony prominences.  Do not use on broken skin - washable/reusable.
	Heel One size	£15.18	2	
	Sheet 10cm x 10cm x 0.3cm	£4.30	1	
	Sheet 10cm x 10cm x 1.2cm	£12.95	1	
	Strips 50cm x 2.5cm x 0.3cm	£5.39	1	
	Strips 30cm x 5cm x 0.3cm	£6.47	1	
Biobag (Larval Therapy)	2.5cm x 4cm (BB50)	£199.50**	n/a	<p><b>For use by specialists in wound care only or under their advice and direction.</b></p> <p>Not for use on bleeding wounds or if close to major blood vessel.</p> <p>Not for use on patients on anticoagulants.</p> <p>Use with caution on wounds over exposed organs.</p> <p><b>**These are 'pharmaceutical specials' so the price is subject to change. Order from manufacturer -</b></p> <p>Telephone: 0845 230 1810 E-mail: orders@biomonde.com Fax: 01656 668 047 <a href="https://biomonde.com/en/">https://biomonde.com/en/</a></p>
	4cm x 5cm (BB100)	£229.50**		
	5cm x 6cm (BB200)	£249.50**		
	6cm x 12cm (BB300)	£279.50**		
	10cm x 10cm (BB400)	£299.50**		

\*Prices from the December 2018 Drug Tariff

## If further information is needed, the tissue viability service for each area can be contacted at:

### Croydon

**Primary care tissue viability team**

Tel: 0208 274 6163

**Hospital tissue viability team** Tel: 0208 401 3000, ext 4204 or bleep 383

### Kingston

**Primary Care, Your Healthcare tissue viability service**

Team email: [tissueviabilityservice@yourhealthcare.org](mailto:tissueviabilityservice@yourhealthcare.org)

Tel via Single Point of Access: 020 8255 1300

**Hospital tissue viability service** email: [khft.tissueviabilityteam@nhs.net](mailto:khft.tissueviabilityteam@nhs.net) or Tel: 020 8546 7711, bleep 135

### Merton

Central London Community Healthcare NHS Trust (CLCH)

Tel: 0333 241 4242

Tel via single point of access: 0333 004 7555

### Richmond

Hounslow and Richmond Community Healthcare (HRCH)

Team email: [HRCH.TissueViability@nhs.net](mailto:HRCH.TissueViability@nhs.net)

Tel via single point of access: 020 8973 3450 or for advice 020 8714 4085

### Sutton

Sutton Health & Care

Tel via single point of access: 020 8661 3908

### Wandsworth

Central London Community Healthcare NHS Trust (CLCH)

Tel via single point of access: 0333 300 2350

**This wound formulary has been developed by the south west London wound-care sub group.**

